# gridinlux

## TRAINER-LEGS REVOLUTION

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IMPORTADOR: B85960425

Dear client,

We would like to thank you for the trust you have placed in our brand, choosing us as a sure guarantee to improve your quality of life.

At Gridinlux we are constantly working to design new products that help to make your day-to-day life easier. We invest an enormous effort into research and development, being able to produce, thanks to the latest technology, efficient systems which are based on the thing that we value the most: people.

From today, we will be with you to help you in every way you need. This is our promise to you.

Remember:

1. The photos of the product in the manual and on the box correspond to the most recent model and only serve as a reference.
2. This manual contains recommendations for use.

**Antonio Castellano**  
European Product Manager

#### Correct disposal of this product

This symbol indicates that, within the EU, this product should not be disposed of along with other household waste. To prevent potential damage to the environment or human health from the uncontrolled elimination of waste, please recycle it responsibly to promote the sustainable reuse of its material resources.

# INDEX

- 22** | Introduction
- 22** | How to use the leg exerciser
- 23** | Program information
- 25** | Remote control functions
- 25** | Usage recommendations
- 26** | Cleaning and maintenance
- 27** | Warranty and Official Technical Service
- 27** | Technical data


## Do you spend hour after hour, day after day sitting or standing without doing any kind of exercise?

This lack of movement may afflict your daily wellbeing, disavouring your positive attitude and feeling of fulfillment.

If you spend long hours each day sitting in the same chair at home or at work, you will have noticed that your legs are weakened and feel uncomfortable after having sat in the same position for a prolonged time.

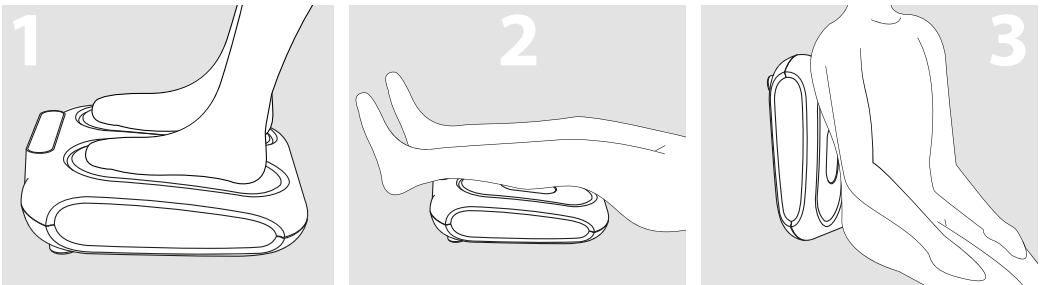The TRAINER-LEGS REVOLUTION foot and leg massager is a motorized solution with the VIBRO SENSE System based on acupressure that will exert a constant and cyclical massage at key points of your body to help you combat that feeling of heaviness and fatigue.

Receive relaxing massages that alleviate the sensation of muscular load in your feet, your legs and even your lumbar:

### How to use the leg exerciser.

- Connect the power adapter to the appliance and plug it into an outlet.
- Sit down and position your feet on the pedals. The default programs The

 POWER	Tap this button to turn on the device and put the machine in standby mode. The standard time it will display is 15 min.
 PROGRAM	With the device on in standby mode, you can select programs HO, UI, P1, P2 and P3.
 START / STOP	Touch this button to start the exercise in MANUAL HO mode or to start the exercise if you have selected a program. Tap this button to stop exercising.
 SPEED  	In MANUAL HO mode, touch these buttons to adjust the speed level from 1 to 10.
 TIME  	Tap these buttons to select the duration time. The maximum setting is 30 min.

### Program information

MODE HO	
VIBRATION AND MASSAGE	CONFIGURABLE TIME AND INTENSITY
MODE UI	
VIBRATION	CONFIGURABLE TIME AND INTENSITY

**P1 - PROGRAM 1**

VIBRATION AND MASSAGE	EVERY MIN SPEED CHANGES
<b>1 - 2 - 3 - 5 - 6 - 7 - 8 - 5 - 3 - 1 - 3 - 5 - 8 - 3</b>	

**P2 - PROGRAM 2**

VIBRATION AND MASSAGE	1 MIN - LEVEL 1
VIBRATION	1 MIN - LEVEL 3
VIBRATION AND MASSAGE	1 MIN - LEVEL 5
VIBRATION	1 MIN - LEVEL 7
VIBRATION AND MASSAGE	1 MIN - LEVEL 9
VIBRATION	1 MIN - LEVEL 10
VIBRATION AND MASSAGE	1 MIN - LEVEL 8
VIBRATION	1 MIN - LEVEL 6
VIBRATION AND MASSAGE	1 MIN - LEVEL 4
VIBRATION	1 MIN - LEVEL 2
VIBRATION AND MASSAGE	1 MIN - LEVEL 1
VIBRATION	1 MIN - LEVEL 4
VIBRATION AND MASSAGE	1 MIN - LEVEL 8
VIBRATION	1 MIN - LEVEL 6
VIBRATION AND MASSAGE	1 MIN - LEVEL 4


### P3 - PROGRAM 3

VIBRATION	EVERY MIN SPEED CHANGES
3 - 5 - 7 - 9 - 2 - 4 - 6 - 8 - 10 - 1 - 5 - 10 - 5 - 1 - 5	

default programs P1-P2-P3 vary the speed and time level of the device automatically.

During the P1-P2-P3 programs, neither the time nor the intensity can be customized. With the UI and HO programs, you can configure both the time and the speed.

### Remote control functions
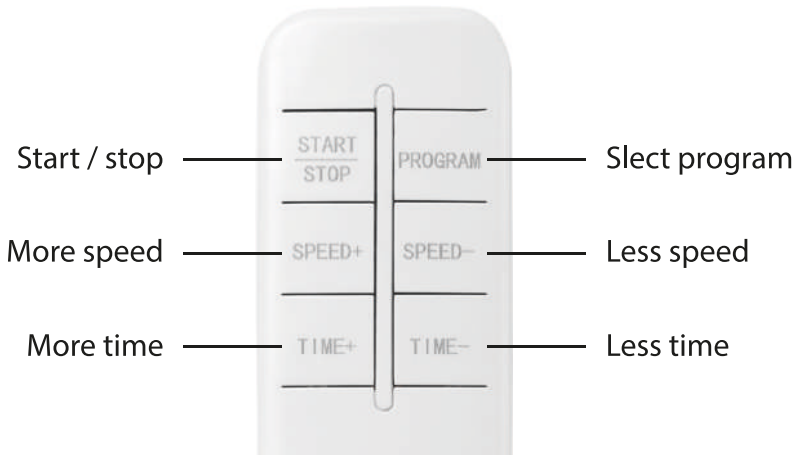
### Usage recommendations.

- Do not sit or stand on the device.
- Keep the device out of reach of children.


- This device does not have any playful purpose and must not be used as a toy.
- Do not use the device on a slippery or inclined surface.
- Do not wet the device.
- Do not put fingers or toes between the pedals.

### Cleaning and maintenance

Before performing any type of cleaning on the device, disconnect it from the network.

You can carefully remove the silicone pads that cover the bases where you rest your feet. So you can clean them without problem. Then put them back carefully, making sure they are 100% dry.

The rest of the device can be cleaned with a damp cloth and a mild detergent. Do not use abrasive or corrosive cleaning agents. Avoid at all times that the device comes into contact with water or other liquids.

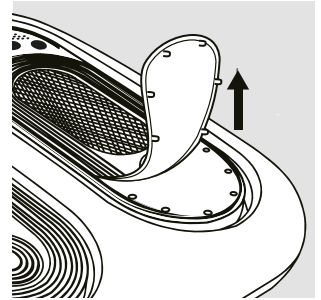TRAINER-LEGS REVOLUTION has been especially designed to improve the wellbeing of people who move very little in their daily lives and have to remain seated for a long time each day.

A constant use in moments when we have to stay seated while performing a parallel task may be the best occasion to use TRAINER-LEGS REVOLUTION and benefit from its positive effect on our sensation of wellbeing.

## Warranty and Official Technical Service

This product is guaranteed for 2 years from the date of purchase, as long as the purchase invoice is well kept and sent to us when needed, the product is in perfect physical condition and is given proper use in accordance with the guidance in this instruction manual .

The warranty will not cover:

- If the product has been used beyond its capacity or usefulness, mistreated, beaten, exposed to humidity, submerged in any liquid or corrosive substance, as well as any other fault attributable to the consumer.
- If the product has been disassembled, modified or repaired by persons not authorized by the official Gridinlux Technical Service.
- If the incident has been caused by normal wear and tear of the parts due to its use over the time.

The warranty service covers all manufacturing defects for 2 years based on current legislation, except consumable parts. In case of misuse by the user, the warranty service will not be responsible for the repair. If on any occasion you detect an incident with the product or have any questions, contact the Official Gridinlux Technical Assistance Service through the email address: [teayudamos@gridinlux.es](mailto:teayudamos@gridinlux.es)

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