# gridinlux

## TRAINER-LEGS SMART

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Dear client,

We would like to thank you for the trust you have placed in our brand, choosing us as a sure guarantee to improve your quality of life.

At Gridinlux we are constantly working to design new products that help to make your day-to-day life easier. We invest an enormous effort into research and development, being able to produce, thanks to the latest technology, efficient systems which are based on the thing that we value the most: people.

From today, we will be with you to help you in every way you need. This is our promise to you.

Remember:

1. The photos of the product in the manual and on the box correspond to the most recent model and only serve as a reference.
2. This manual contains recommendations for use.

**Antonio Castellano**  
European Product Manager

#### Correct disposal of this product

This symbol indicates that, within the EU, this product should not be disposed of along with other household waste. To prevent potential damage to the environment or human health from the uncontrolled elimination of waste, please recycle it responsibly to promote the sustainable reuse of its material resources.

# INDEX

- 22** | Introduction
- 23** | How to use the leg exerciser
- 24** | Program information
- 24** | Remote control functions
- 25** | Usage recommendations
- 25** | Exercise on a daily
- 27** | Warranty and Official Technical Service
- 27** | Technical data

## **Do you spend hour after hour, day after day sitting or standing without doing any kind of exercise?**

This lack of movement may afflict your daily wellbeing, disavouring your positive attitude and feeling of fulfillment.

If you spend long hours each day sitting in the same chair at home or at work, you will have noticed that your legs are weakened and feel uncomfortable after having sat in the same position for a prolonged time.

The personal leg exerciser TRAINER-LEGS SMART is the only motorized solution of a passive exerciser that creates a constant leg movement to stimulate your lower extremities while you remain seated. You can use this device with or without shoes, to make the most of its functions.

It allows you to comfortably train your legs, while traditional exercise techniques require space and positions that, due to various difficulties, for many people are impossible to realise. With TRAINER-LEGS SMART you also benefit from the properties of magnetism.

## How to use the leg exerciser.

- Connect the power adapter to the appliance and plug it into an outlet.
- Sit down and position your feet on the pedals.

 POWER	Tap this button to turn on the device and put the machine in standby mode. The standard time it will display is 15 min.
 PROGRAM	With the device on in standby mode, you can select programs P1, P2 and P3. In these programs the level varies automatically, and you can select the time between 1 and 30 min.
 START / STOP	Touch this button to start the exercise in MANUAL HR mode or to start the exercise if you have selected a program. Tap this button to stop exercising.
 SPEED 	In MANUAL HR mode, touch these buttons to adjust the speed level from 1 to 30.
 TIME 	Tap these buttons to select the duration time. The maximum setting is 30 min.

Smooth and cyclical movements.

If you want to tilt the appliance, remove the legs from the bottom.

**Program information**

<b>HR</b>	Manual mode: configurable in time and speed
-----------	---

In MANUAL HR mode, the time can be set between 1 and 30 min, and levels from 1 to 30 can be chosen.


<b>P1</b>	<b>1 2 3 4 5 6 7 8 9 10 9 8 7 6 5</b>
<b>P2</b>	<b>11 12 13 14 15 16 17 18 19 20 19 18 17 16 15</b>
<b>P3</b>	<b>21 22 23 24 25 26 27 28 29 30 29 28 27 26 25</b>

The default programs P1-P2-P3 automatically vary the speed level of the device.

You can choose the duration of the programs between 15 and 30 minutes. The programs work as follows: increase 10 speed levels in 10 minutes and decrease 5 speed levels in 5 minutes.

If the selected run time is longer, the cycle starts again.

**Remote control functions**
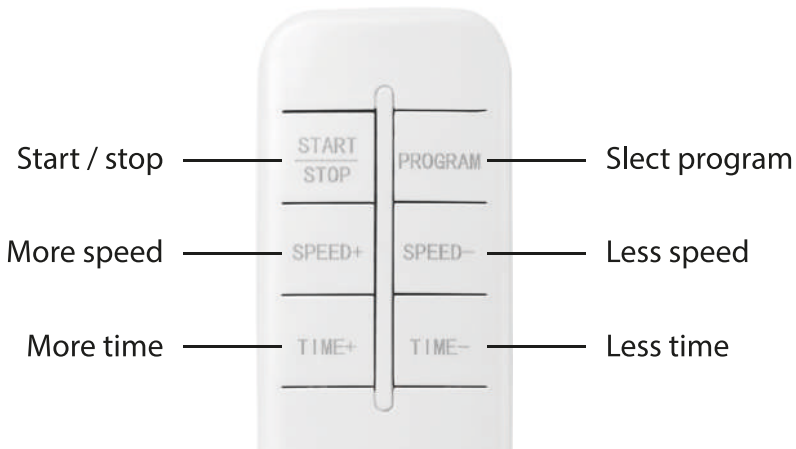
## **Usage recommendations.**

- Do not sit or stand on the device.
- Keep the device out of reach of children.
- This device does not have any playful purpose and must not be used as a toy.
- Do not use the device on a slippery or inclined surface.
- Do not wet the device.
- Do not put fingers or toes between the pedals.

## **Exercise on a daily basis to improve your wellbeing.**

TRAINER-LEGS SMART has been especially designed to improve the well-being of people who move very little in their daily lives and have to remain seated for a long time each day.

A constant use in moments when we have to stay seated while performing a parallel task may be the best occasion to use TRAINER-LEGS SMART and benefit from its positive effect on our sensation of wellbeing.

More and more professionals find themselves forced to remain seated during their entire working day. The leg exerciser is the best tool to prevent discomfort and pain caused by a lack of exercise.

Thanks to this device, your legs stay active because of the constant leg movement it creates. This movement visibly improves the positive feeling of our body.

In addition to that, since device is a passive exercising system, you will train your legs without any effort.

For example, 3 hours of exercise have a similar result as 1 hour of walking.

## Warranty and Official Technical Service

This product is guaranteed for 3 years from the date of purchase, as long as the purchase invoice is well kept and sent to us when needed, the product is in perfect physical condition and is given proper use in accordance with the guidance in this instruction manual .

The warranty will not cover:

- If the product has been used beyond its capacity or usefulness, mistreated, beaten, exposed to humidity, submerged in any liquid or corrosive substance, as well as any other fault attributable to the consumer.
- If the product has been disassembled, modified or repaired by persons not authorized by the official Gridinlux Technical Service.
- If the incident has been caused by normal wear and tear of the parts due to its use over the time.

The warranty service covers all manufacturing defects for 3 years based on current legislation, except consumable parts. In case of misuse by the user, the warranty service will not be responsible for the repair. If on any occasion you detect an incident with the product or have any questions, contact the Official Gridinlux Technical Assistance Service through the email address: [teayudamos@gridinlux.es](mailto:teayudamos@gridinlux.es)

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